


**Protectors for hip stem cone  
STANDARD Datasheet**


These pieces were designed to protect the cones during production and during installation. Cones protector can be removed after the introduction of hips on the patient avoids the risk of damage during surgery.

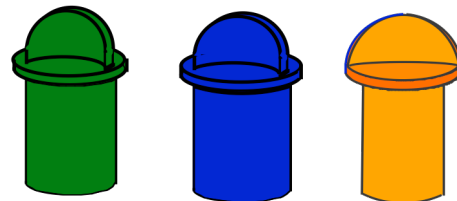

**STANDARD sizes : 10/12 - 11/13 - 12/14**

Raw material is SANTOPRENE ref- 8281-75 MED . medical quality.

Production by GRANDUPLA Fully qualified ISO 14644-1 :1999(E)  
Production controlled atmosphere  
Pack in double PE bag.

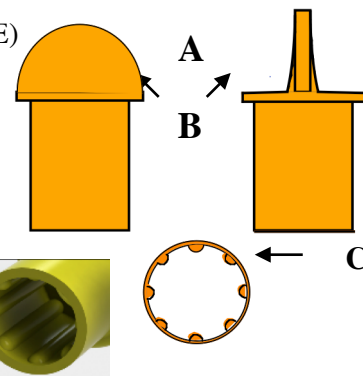

**Two possible grip:**

- A – The upper part with a tool
- B – Part overflowing for a handgrip

Inside Gadroons (C) limits the friction with the stalk of hip. It is also studied to facilitate the introduction or remove the protector after the operation of pose on the patient..

**The main advantages:**

- **Protection of the cone of the stalk of hip during the finish**
- **Protection of the cone of the stalk of hip during the pose, facilitating the work of the surgeon**

We can make special parts, by creating only cavities adapted to our equipment of injection, if you wish a shape or a specific adaptation.

